# Molift Raiser Pro

***molift***<sup>®</sup>  
by Etac

EN - User manual

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SV - Bruksanvisning

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NO - Brukermanual

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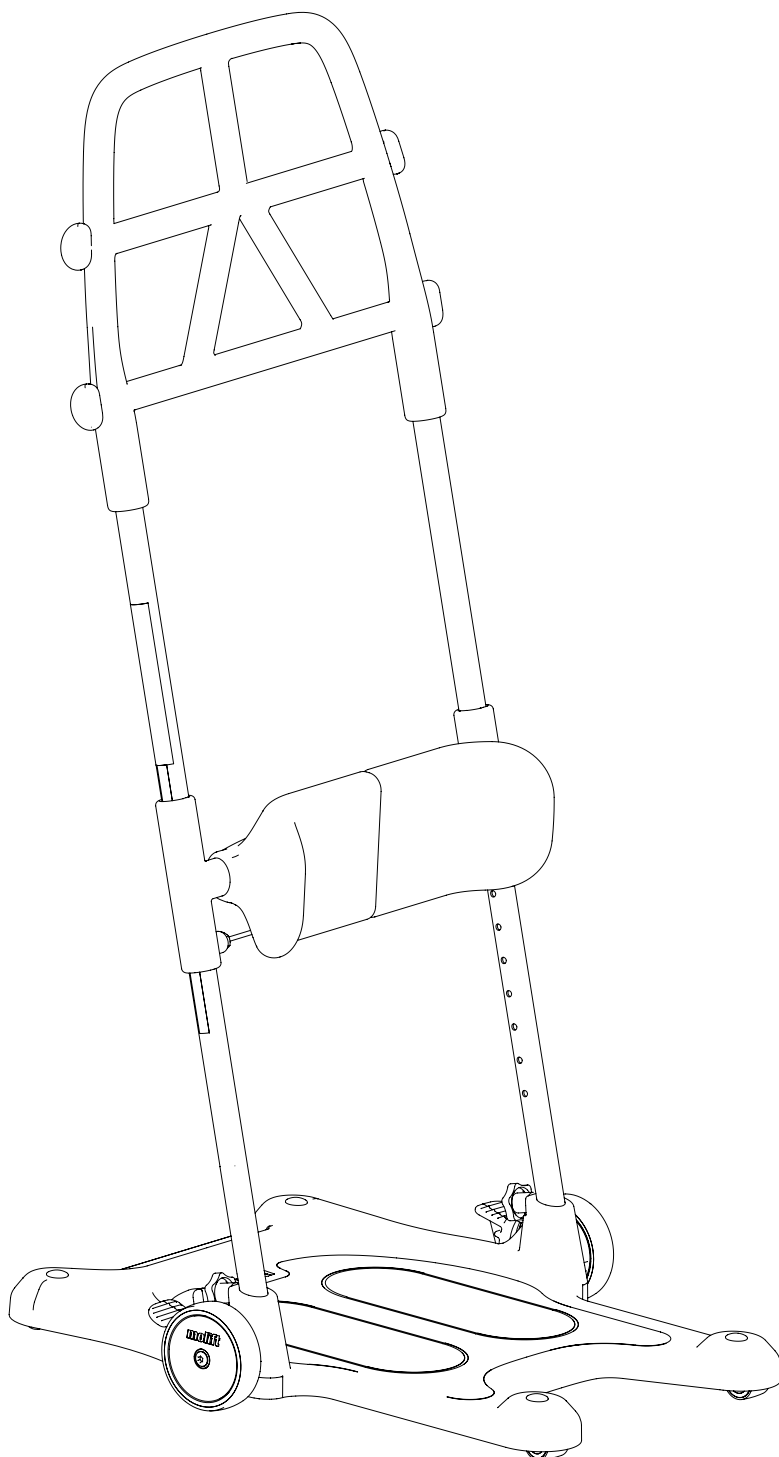
BM43099 Rev. I 2024-01-12

## Contents

User manual - English	- 4
Bruksanvisning - Svenska	- 9
Brukermanual - Norsk	- 14
Brugsvejledning - Dansk	- 19
Käyttöohje - Suomi	- 24
Gebrauchsanweisung - Deutsch	- 29
Handleiding - Nederlands	- 34
Manual d'utilisation - Français	- 39
Manual utente - Italiano	- 44
Manual de usuario - Español	- 49

## General

### CE declaration

The product, and its related accessories, described in this instruction for use (IFU), is in compliance with the regulation (EU) 2017/745 of 5. April 2017 – as a medical device, risk class I.

The product is manufactured and tested according to EN 12182:2012 and selected parts of EN/ISO 10535:2006.

Any serious incident that occurs in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

### Conditions for Use:

Transferring a person may pose an element of risk. Only trained personnel should use the equipment covered by this user manual.

### Warranty notice:

Two-year warranty against defects in workmanship and materials of our products. Please refer to [www.etac.com](http://www.etac.com) for terms and conditions.

Visit [www.etac.com](http://www.etac.com) to download the latest version of our documentation.

## Important

Read User Manual before use!

It is important to fully understand the content of the user manual before attempting to use the equipment. This User Manual contains important safety instructions and information regarding the use of the Molift Raiser Pro. For instructions/information regarding the accessories, see separate user manual.

In this manual the user is the person being transferred.

The assistant is the person operating the Molift Raiser Pro.

### Modifications and use of components made by other manufacturers:

We recommend only using Etac components and spare parts. Declaration of conformity is not valid and Etac is not responsible for warranty if any modifications are made to the Product. Etac shall not be liable for faults or accidents that can occur when using components made by other manufacturers.

## Risk assessment

It is important that the Molift Raiser Pro has been tested with the individual user and for the intended transfer. Assess the risks and take notes. Evaluate the user's ability to use the Product and if needed use any suitable Molift Raiser Pro accessories. You as a carer are responsible for the safety of the user.

Only technical personnel are allowed to perform repairs.

## Introduction

The person referred to as "the user" is the person who is standing on the product. "The carer" is the person who manoeuvres the product.

## Intended Purpose

The sit-to-stand aid is an assistive device intended for alleviation of, or compensation for, a functional impairment due to an injury or disability. The device is designed for an individual lacking the ability to stand up and transfer themselves over shorter distances to another sitting position to/from a bed, a wheelchair, a chair, a toilet or similar due to reduced mobility or physical strength.

## Intended User

The device can be used for children, adults and elderly people with a need for support to stand up and be transferred over shorter distances between sitting positions.

## Intended Environment

For indoor use on horizontal surfaces in acute care, long-term care and home care.

## Indications

For users who are in need of support to stand up and be transferred over shorter distances between sitting positions; e.g. people suffering from pain, reduced range of motion in their joints, stiffness and muscle weakness. Not an exhaustive list.

The user group for the device is based on individual health and mobility function and not on a specific diagnosis or age.

## Contraindications

There are no known contraindications.

## Precautions

The user should be able to follow instructions, grip and pull up their own body weight, maintain a good core stability and load at least one leg.

## Training Requirements

Only carers who have read and understood the user manual are allowed to use the device.

Labels and symbols

Product label:

The label contains art. number and serial number (bar codes).  
The EAN code contains;

(01) (EAN no.): xxxxxxxxxxxxxxxx

(11) (Batch): YYMMDD

(21) (Serial number): xxxxx

Symbols:

- Manufacturer
- Date of manufacture  
YYYY-MM-DD
- Catalogue number
- Serial number
- Medical Device
- CE marked
- Refer to user manual
- Max temperature/time when washing.  
max 70°C  
10 min
- Max user weight (SWL)  
max 170 kg  
max 375 lbs
- pH when washing  
min 5  
max 9
- Important safety information
- QR code

Short user guide:

The short user guide graphic includes a QR code, an information icon, a wheelchair icon, and a diagram of the device. Below the diagram is the text:  
**Raiser Pro**  
**molift**  
 by Etac  
 max 170 kg  
 max 375 lbs  
 max 70°C  
 10 min

Technical data

Safe Working Load (SWL):

170kg (375 lbs)

Weight of unit:

10kg

Dimensions:

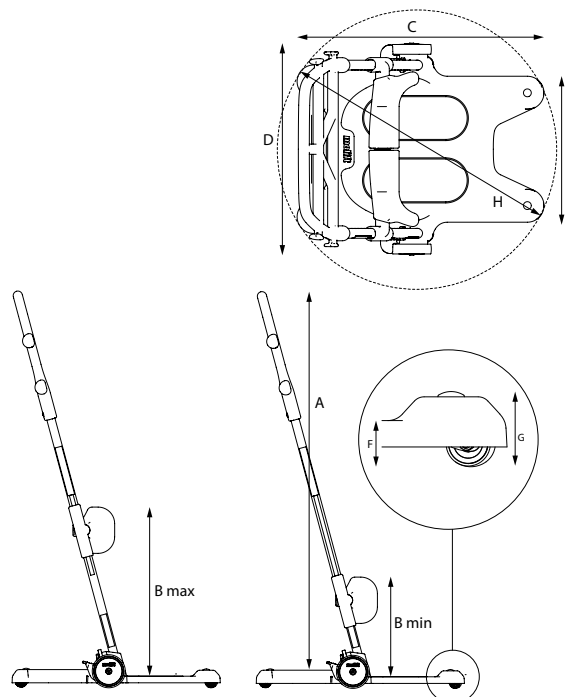
- Height(A): 1,2 m
- Height (B): 0,3-0,51 m
- Depth (C): 0,63 m
- Width (D): 0,54 m
- Width (E): 0,37m
- Step in height (F): 0,035m
- Clearance (G): 0,057m

Material:


Powder coated aluminum,  
Plastic, Steel.

Expected Lifetime:


The Molift Raiser Pro has an expected lifetime of 5 years or 15000 transport cycles.

Supplied parts

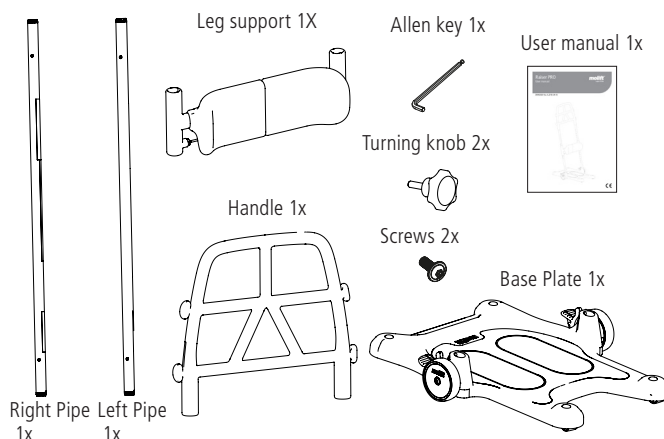
Molift Raiser Pro

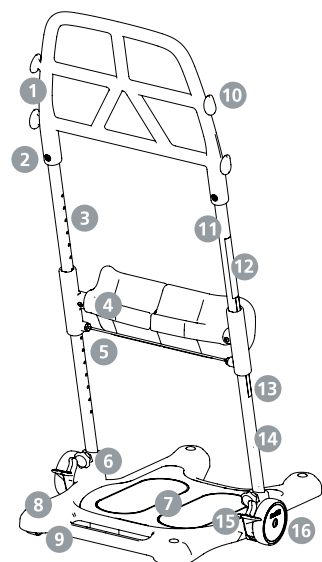
1. Handle
2. Screw
3. Left Pipe
4. Leg support
5. Release wire
6. Knob
7. Foot plate
8. Base plate
9. Castors Ø35
10. Attachment for accessories (2 heights)
11. Right Pipe
12. Short user guide
13. Height label
14. Product label
15. Brake pedal
16. Wheel Ø100

### Mounting pipes to base plate

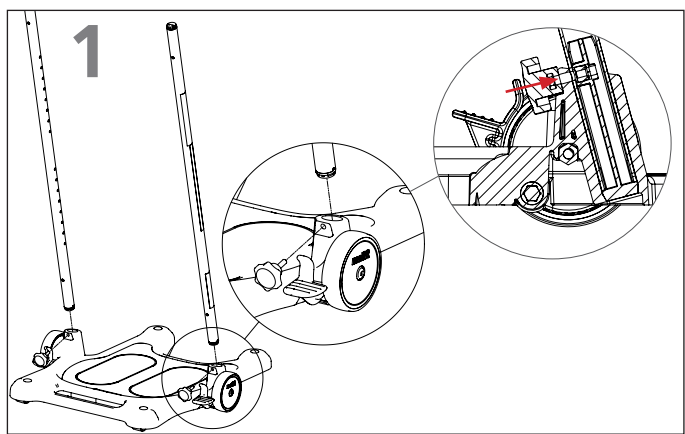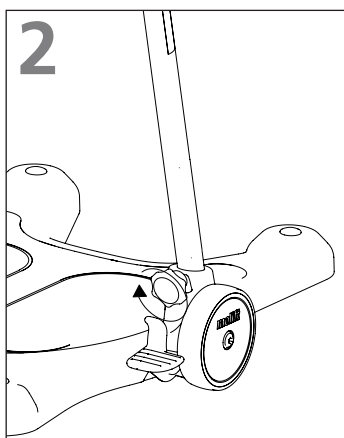1. Place the base plate on an even surface and lock brakes. Insert pipe with holes facing towards center of plate. Pipe with label must be at the right, standing on the assistant side. Mount the knob so it is perpendicular to the pipe and angled slightly up in relation to the horizontal.


2. Tighten knob firmly by turning knob clockwise.


**Make sure that knob are Properly fastened. Loose knobs can make the pipes wobbly and the Molift Raiser Pro unsafe to use.**
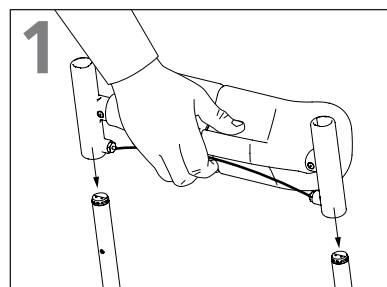
### Mounting leg support

The knee cushions are adjusted individually. The leg support has 7 step adjustment from 0,3 to 0,51 meters. Follow these two steps when mounting the leg support:

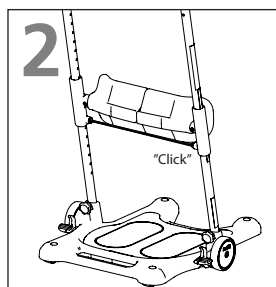
1. Pull release wire and slide the leg support on pipes, cushion facing user side.


2. Release wire and let the leg support settle in desired position.

### Mounting handle

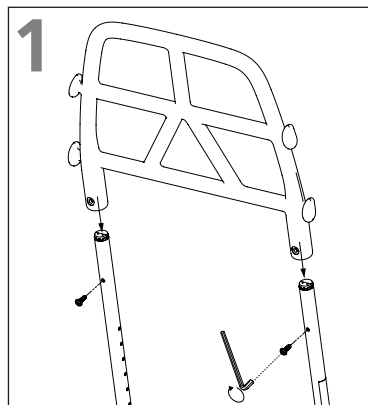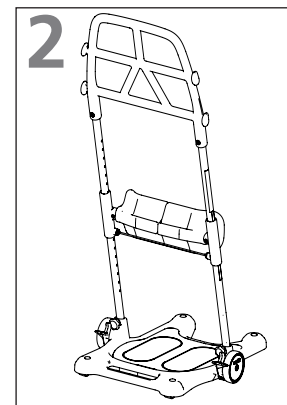1. Slide the handle on pipe so holes on pipe matching holes on handle. Insert unbraco screws and tighten with allen key.

2. Make sure unbraco screws are pointing towards assistant side.

### Checklist after assembly

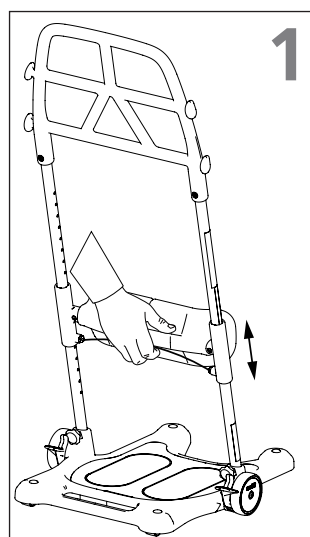
**Don't start using Molift Raiser Pro before completing the checklist after assembly.**


Use this checklist to verify that the Molift Raiser Pro is Properly installed and can operate correctly and safely before use. See assembly information on the previous page for details.

- Make sure that the leg support is mounted correctly.
- Make sure that both knobs between outer pipes and base plate are mounted correctly and Properly tightened
- Make sure that the handle is mounted correctly.
- Make sure Molift Raiser Pro moves freely.
- Make sure brake works.

### Adjusting leg support height

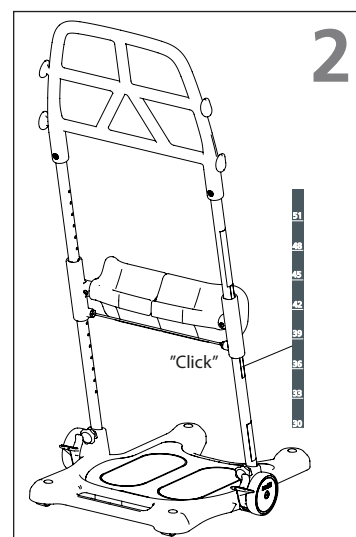1. Pull wire on the back side of cushions. Move leg support up or down to desired height.


2. When correct knee support height is obtained, release wire and let the leg support settle in desired height. You will hear the leg support "click" in place and the red wire is straight. The knee support height can be read from the scale on the pipe. The number above the leg support bracket refers to the correct height.


 Make sure both sides are properly settled. Push up or down until you get the auditory feedback click. Unattached leg support can cause knee support to slide down and make the Molift Raiser Pro unsafe to use.

 The wheels must always be locked when the user stands up and sits down.

**Positioning Molift Raiser Pro**

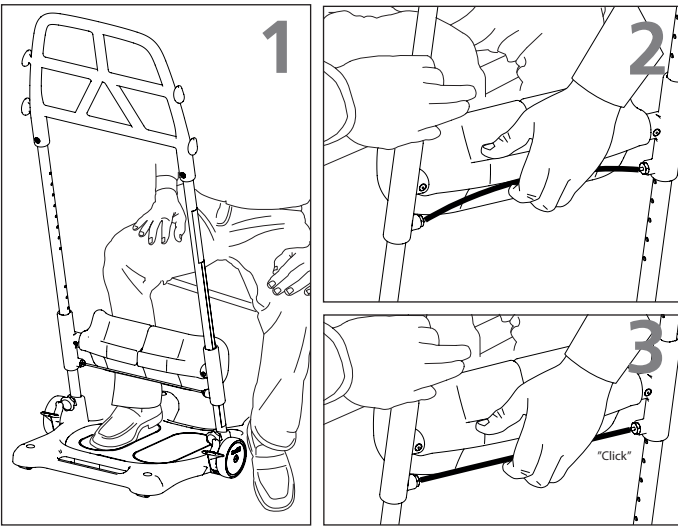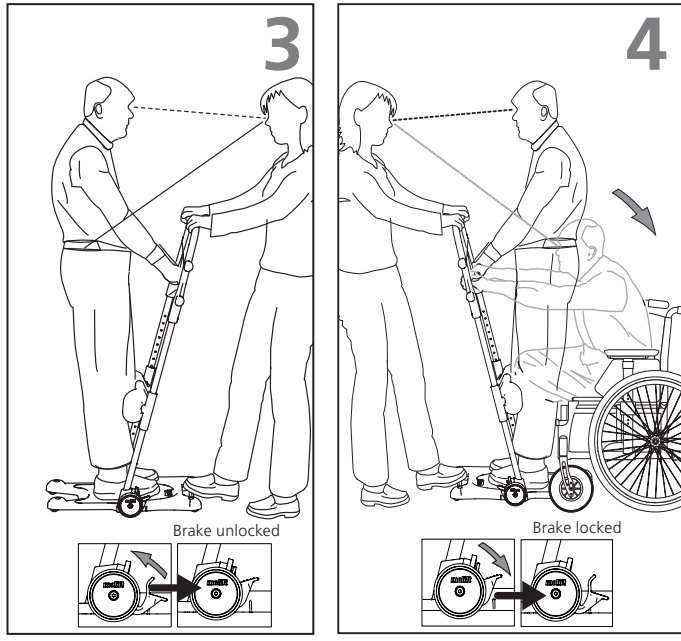
1. The user's feet are placed on the footplate.
2. Adjust the height of the leg support so that the top of the support is below the knee of the user.
3. When desired height is reached the leg support "click" in place.

3. Release the brake with one foot. Pull the Molift Raiser backwards and then steer it toward the new seating position. Make sure that any wheels are locked in the new seating. The transfer is best done with the user faced in the direction of travel.
4. Lock both wheels. Place a foot on the edge of the footplate before the user sits down. Counter balance while the user sits down. Maintain eye contact with the user during the entire transfer Process.

**Transferring user**

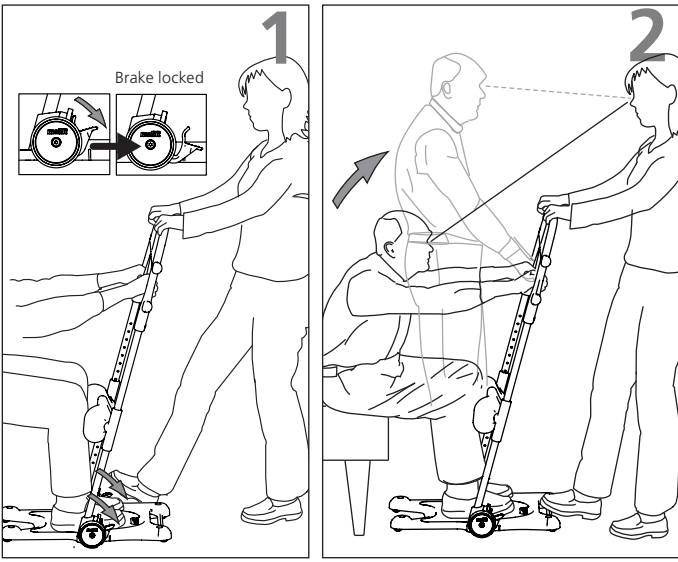
1. Lock both wheels by pressing down the brake on one side. Both wheels lock regardless of which side assistant steps on. Place one foot on the edge of the base plate. Choose a grip on the handle that is comfortable for both assistant and the user.
2. Encourage the user to stand up, at the same time, counter balancing the users weight using a foot and hands (use your body-weight by leaning backwards).

 The carer must always counter balance the weight of the user when she/he stands up or sits down.


**Transport and storage**

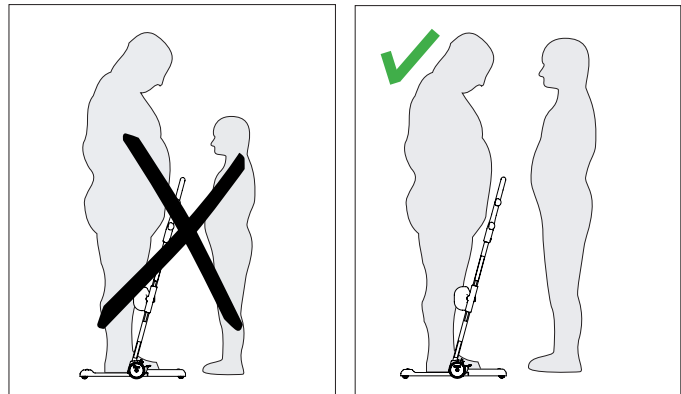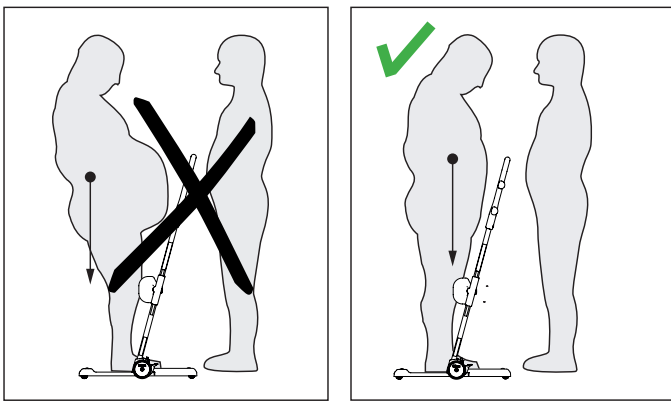Molift Raiser Pro can be stored and transported under temperatures between -25 - 70°C and relative humidity between 15 - 93%. Don't use the Molift Raiser Pro for transportation of goods.

**User / Assistant**

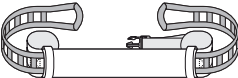
The assistant has to be able to handle the weight of the user.

Make sure the users weight is close to the center of the base plate.  
To avoid serious injury, the user shall avoid extended reaching or leaning movements.

### Accessories

-  Molift Raiser Soft Handels
-  Molift Raiser Strap
-  Molift Raiser Strap +
-  Molift Raiser Belt
-  Molift Raiser Pro Heel Strap
-  Molift Raiser Pro Block

### Cleaning and disinfection

Clean on a regular basis. Clean surfaces with a damp cloth using an appropriate pH-neutral detergent. Do not use solvents or strong liquids, this may damage surfaces on the Molift Raiser Pro. For disinfection when needed; use isoPropyl alcohol or similar. Avoid abrasive cleaning Products. Remove hair and pile from the wheels, verify that the wheels rotates freely and brakes functions Properly. Make sure not to damage or remove labels when cleaning.

### Reconditioning

Follow cleaning and assembly Procedure. Separate reconditioning document is available at Etac.com.

### Recycling

Molift Raiser Pro can be disposed of as metal and plastic waste. Recycling instructions is available at Etac.com

### Spare parts

A list of spare parts is available at Etac.com.

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Etac A/S