







Multi Bari

Designed for comfort and pressure management control, the Multi Bari has removable waterfall back cushions that can be individually positioned and re-positioned for comfort with the option to add or remove fibres to alter their firmness. Couple this with individually positioned back blockers offering further lumbar adjustments, plus the option for the seat and width measurements to be built to the user's requirements, and this chair is ideal for people seeking deep comfort and enhanced posture control.

The Dual Motor mechanism with rise function, and overall 254 kg (40 stone) weight limit, and elevated leg rest weight limit of 160 kg (25 stone), ensures maximum support.

Key features at a glance:

- Standard Dual Motor rise and recline mechanism
- Adjustable and removable waterfall back cushions
- Back blockers for additional lumbar support
- Choice of interchangeable pressure management seat cushions
- Made to measure seat height, width and depth
- Flat wide arms for deep comfort and support
- Seat and leg rest generously padded with memory foam and removable cover
- Loop for handset control and deep side pocket
- Heavy duty glides are built in as standard.

Additional options:

- Full length arm covers
- Head and seat covers
- Replacement chaise cover and memory foam
- Additional adjustment blockers for back cushions
- Enclosure skirt
- Range of headrests.

Featured fabrics:

Outer arms and wings in Highland Navy, inside arms and wings in Lime and inner in blue vapour permeable multi-way stretch fabric.



Practical / Adaptable and Deeply Comfortable

At Repose we pride ourselves on delivering high levels of comfort, posture control and pressure management to support an individual's health and lifestyle options. Our plus size Multi Bari chair is no exception with the following key comfort features:



Waterfall back cushions delivering adjustable comfort

The three removable waterfall cushions can be repositioned enabling the user to be seated and supported posturally to meet their individual needs. The comfort and postural options are further enhanced with the addition of blockers that adjust the shape and support of the back cushions, and the ability to adjust the firmness of each individual cushion by adding or taking away fibres.

This back style can be useful for those with mild kyphosis, scoliosis and reduced muscle mass on their back.

Seat Dimensions – Made to Measure

The Multi Bari comes with a range of made to measure options:

Seat Height ensuring that the client's feet are flat on the floor and that their thighs are in contact with the full length of the seat cushion to provide a comfortable sitting position.

Seat Width supporting correct pelvic positioning and stability as well as reducing any pressure areas especially around the hip joint. Always allow adequate room for natural movement and different changes of clothing.

Seat Depth ensuring that the thighs are supported throughout their length, but not too long or too short so as to avoid pressure points at the back of the knees and on the ischial Tuberosities.


Kate Sheehan
Independent Occupational Therapist



Interchangeable seat cushion and chaise

There is a choice of four pressure management seat cushions all designed to provide different options. This enables the chair to be adapted as the needs of the user change.

In addition, there is a full seat cover and chaise to reduce the risk of pressure areas developing due to pinch points and to reduce shear and friction during the chair movements.

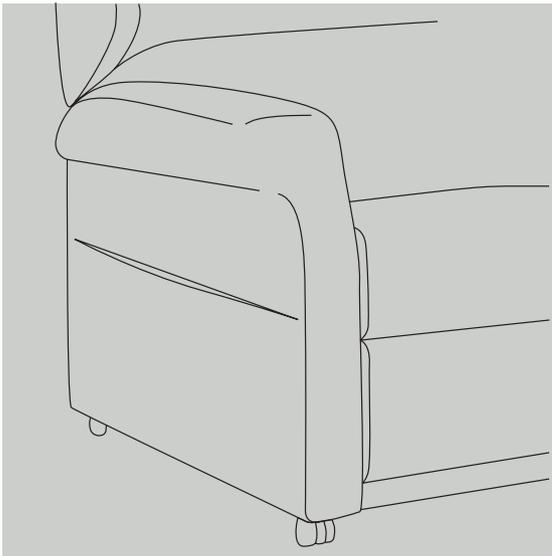

Arm Rest Design

The arms have been designed to be extra flat and wide, this supports upper limb weight and reduces strain on joints and muscles.

Standard Dual Motor Rise and Recline Mechanism

The Multi-Bari comes with a Standard Dual Motor rise and recline mechanism which is ideal for those clients:

- who are unable to independently shift their weight either in a reclined or seated position,
- with poor skin integrity where there is a high risk of skin breakdown.
- with complex needs that require carer support to change position without putting their own health at risk whilst carrying out positional changes.



Full Rise



Static Position



Independent Back Rest



Semi Recline



Footrest Horizontal



Fully Reclined

Standard Dual Motor

The two motors allow the back rest to move independently of the leg rest enabling a wide range of relaxing positions, from sitting upright and raising the leg rest to relieve tired legs, to a fully reclined resting position. Leg section elevation can support reduction in oedema and fluid retention. This reduced fluid in the lower limbs allows a user to maintain their standing and walking abilities.

When the seat is in the rise position it tilts transferring the user's weight to their knees and taking away the pressure on the hips allowing them to walk away from the chair more easily.

More information and clips showing this mechanism in action are available on our website.

Adaptable Back and Seat Cushions

For ease of cleaning and replacement all back and seat cushions and covers are removable.



Waterfall back

The Multi-Bari comes with our deeply comfortable three cushioned, waterfall back, all fibre filled. They can be easily repositioned on the back for the client's individual requirements and can also have their firmness adjusted by adding or taking away fibres.



Back blockers for postural support

At the back of the waterfall cushions, sleeves are fitted to accommodate the blockers to provide additional postural support. The position of each back cushion can be adapted with blockers to adjust the shape and support for the specific needs of each user. Blockers are provided as standard.



Seat cushion and chaise removal


Seat cushions can be removed and interchanged to adapt the chair to each new user's pressure care or support needs.

Easy to change covers. Quick to remove and replace as demonstrated on the removable chaise.

Pressure Management Seat Cushions

The Multi-Bari comes with our interchangeable seat cushions. Purchase one or more cushions to support the different needs of the individuals who may use this chair. With all these seat cushions, we recommend a fabric with four way stretch properties to aid the reduction of shear and friction.

Viscoform

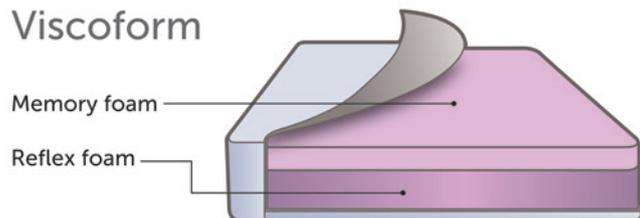

Memory foam


Reflex foam

Pressure Risk

LOW MEDIUM

For those who have good posture and mobility but are seated for greater than one hour at a time and require pressure relief and supported comfort. This cushion has a high-density foam base which offers a supportive and firm seat base with a top layer of memory foam to allow immersion into the cushion, giving stability and a large surface area to reduce pressure levels.

Coolform

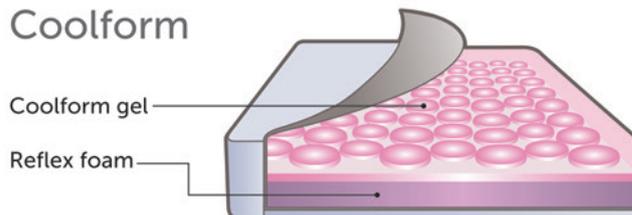

Coolform gel

Reflex foam


Pressure Risk

MEDIUM HIGH

For those sitting for long periods with specific comfort issues such as a high risk of pressure sores.

This cushion uses the latest medical grade silicone gel technology allowing pressure to be evenly distributed. The gel conforms to the shape of the body and equalises to body temperature. The gel cells allow air to flow easily. Coolform is totally odourless and hypoallergenic.

Liquiform

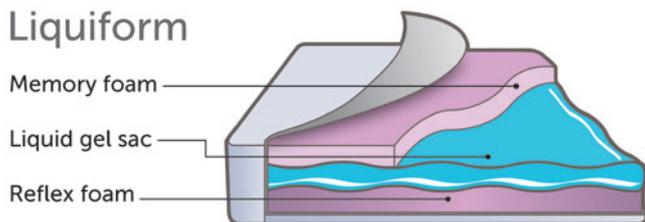

Memory foam

Liquid gel sac

Reflex foam


Pressure Risk

MEDIUM HIGH

For those who have good posture and mobility but sit for long periods of time requiring pressure relief and supported comfort. Recommended for people with limited muscle bulk or at risk of excessive weight loss.

This cushion has a high-density foam which offers a supportive and firm seat base. The foam is topped with a fluid liquid gel in a soft silicone cover which provides a high level of pressure relief. The liquid gel helps to dissipate heat and constantly adapts to the user's movement.

Dynaform*

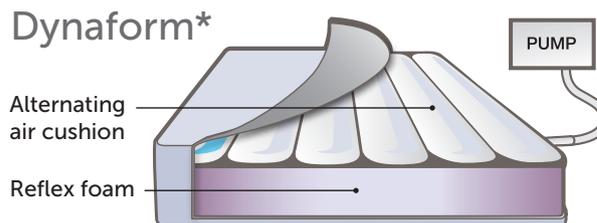

Alternating air cushion

Reflex foam

Pressure Risk

VERY HIGH

For those at very high risk of developing pressure sores or a history of sores where the highest priority is to prevent further pressure damage.


The Dynaform cushion is a dynamic alternating air cushion, which constantly changes the air levels within its air cells to recreate the effect of normal body movements, which provides pressure redistribution to prevent pressure build up.

Optional Accessories


Headrests

A range of headrests are available with the Multi-Bari:


- Head roll
- Small head pillow
- Small profile headrest
- Large profile headrest
- Midline headrest.

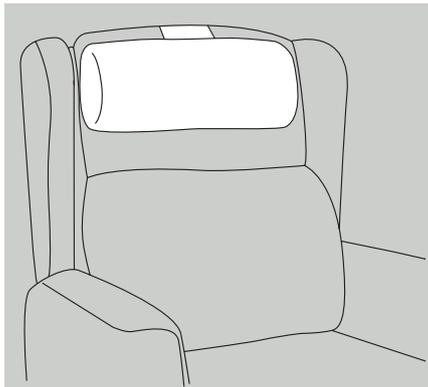


Head roll

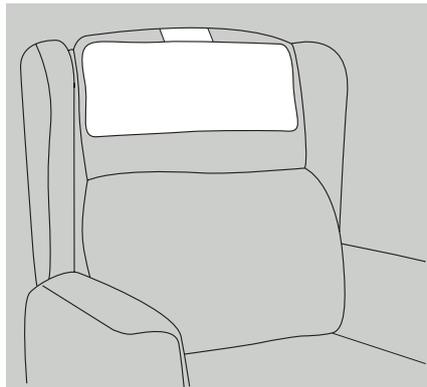


Small head pillow

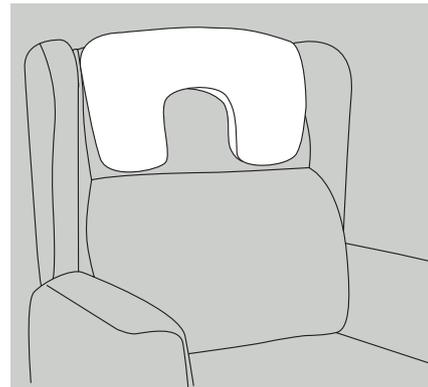

Small profile headrest

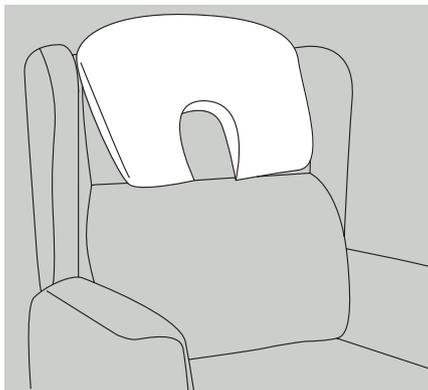


Large profile headrest

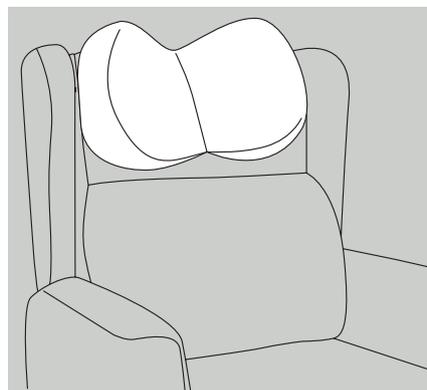


Midline headrest

Optional Accessories (continued)


- Full length arm covers
- Head and seat covers
- Replacement chaise cover and memory foam
- Back blockers
- Enclosure skirt.

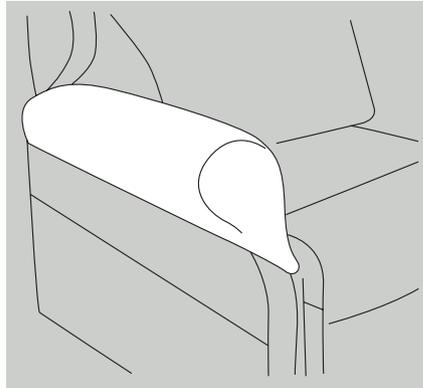

Full length arm cover

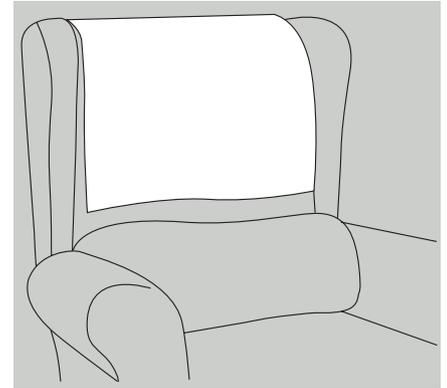

Head cover

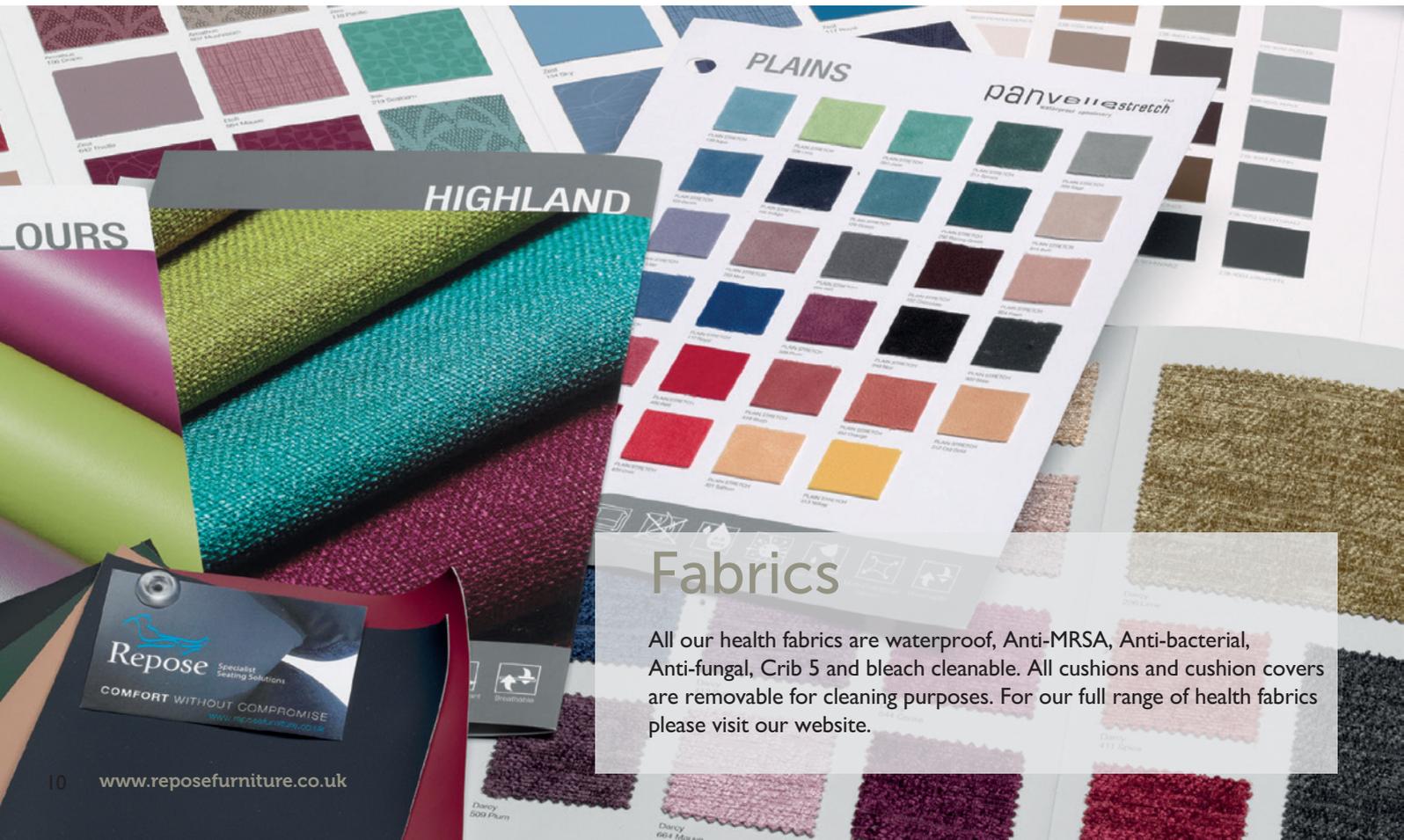

Fabrics

All our health fabrics are waterproof, Anti-MRSA, Anti-bacterial, Anti-fungal, Crib 5 and bleach cleanable. All cushions and cushion covers are removable for cleaning purposes. For our full range of health fabrics please visit our website.

Weights and Dimensions

The Multi-Bari has a weight limit of 254 Kg (40 stone) and a leg rest with a capacity of 160 Kg (25 stone) as standard.



Seat Height	457mm (18")	457mm (18")	457mm (18")	457mm (18")
	508mm (20")	508mm (20")	508mm (20")	508mm (20")
Seat Width	610mm (24")	660mm (26")	711mm (28")	762mm (30")
Seat Depth	457mm (18")	457mm (18")	457mm (18")	457mm (18")
	508mm (20")	508mm (20")	508mm (20")	508mm (20")
	559mm (22")	559mm (22")	559mm (22")	559mm (22")

Please note that depending on the seat cushion built into the chair, the seat height may vary slightly from that shown above.

Warranty

Two years on mechanisms, frame and upholstery.

One year on air seat cushions.

Full terms and conditions are available on our website.

We recommend that professional advice is sought when purchasing specialist seating.



Repose Furniture Limited
Attwood House
Cokeland Place
Cradley Heath
West Midlands
B64 6AN

Tel: 0844 7766 001*

Fax: 0844 7766 002*

Email: info@reposefurniture.co.uk

*Calls cost 7p a minute plus your networks access charge.



www.reposefurniture.co.uk



Multi Bari